


User Guide

CD81 ExoLISA

Exosome Detection System

Europium time-resolved immuno-
fluorescence assay for detection
of exosome antigens

Cat EX503

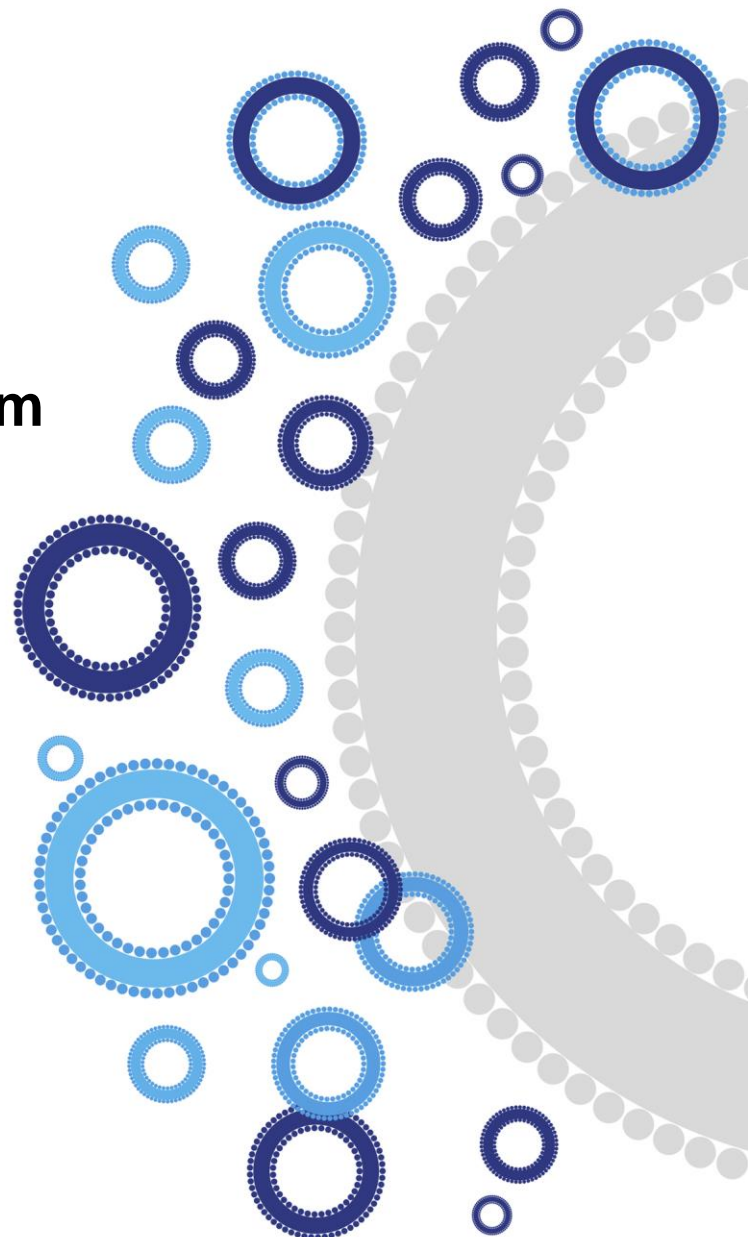

Contents

Storage	4
Product components	4
Equipment and materials required but not supplied with this kit.....	4
Product information	5
Introduction and assay principle	5
General note.....	6
Reagent preparation	6
Technical positive control preparation.....	6
See example assay reading on page 9.	7
Protocol.....	7
Coat the wells with CD81 capture antibody	7
Add the sample	8
Add the europium-labeled CD81 detection antibody	8
Signal enhancement and reading	8
Technical positive control example	9
Related products.....	9
Exo-spin™ purification kit.....	9
NTA size profiling service	9
Purchaser Notification.....	11

CD81 ExoLISA

Exosome Detection System

Storage

- Store all components at 4°C.
- The kit has a shelf life of at least three months from receipt.

Product components

- 1 x Streptavidin-coated 96-well plate (8-well strip format)
- 1 x Technical positive control, lyophilized, 25 µg.
- 1 x Eu-labeled CD81 mAb (in 40 µl TSA buffer, 0.1% BSA), 2.75 µg
- 1 x Biotinylated CD81 mAb (in 22 µl PBS 7.4 pH, 15 mM NaN₃), 22 µg

The CD81 antibody used in the kit is human specific.

- 1 x Assay buffer, 22 ml
- 1 x Europium fluorescence intensifier (EFI Solution), 11 ml
- 1 x 25x wash buffer, 20 ml
- 1 x User Guide

Equipment and materials required but not supplied with this kit

- Time-resolved fluorescence microplate reader
- Automatic plate washer
- Plate shaker
- Pipettes for dispensing reagents
- Multichannel micropipette reservoir
- Distilled / Milli-Q water
- Phosphate buffered saline (PBS)

Product information

Introduction and assay principle

In the ExoLISA exosome assay, the same antibody is used for binding of target to the assay plate and for detection. The assay consists of a monoclonal antibody (labelled with biotin) bound to a streptavidin-coated plate which captures protein present on the surface of exosomes (Figure 1). Subsequently, an identical monoclonal antibody (labelled with europium) is used for detection. Because the capture and detection antibody are identical, they require two linked copies of the same epitope for a signal to be detected. Exosomes provide an ideal structure to link CD81 molecules and allow detection of CD81 in this assay. Exosomes typically have multiple copies of CD81 facing towards the attachment surface and additional CD81 molecules available for detection. Any non-specific binding of capture and detection antibodies is unlikely to generate a signal. Using a europium fluorophore (see below) provides high levels of sensitivity for the assay, which is able to detect small changes in the abundance of the target CD81 protein even within unpurified complex biological samples, such as blood plasma and cerebral spinal fluid.

Fluorophores are chemical substances that emit light following excitation by light or other electromagnetic radiation. The emission of light from a fluorophore is maximal immediately following excitation and decays over a period of time. Time-resolved fluorimetry uses fluorophores which have long decay periods. For such fluorophores, measurement of emitted light can be performed when the excitation light is no longer present, thus increasing sensitivity.

Europium is a fluorophore which produces an extended emission decay and has a wide Stokes shift with maximal excitation at 340 nm and peak emission at 615 nm. ExoLISA assays are time-resolved immunofluorescence assays which utilize europium and have been developed to measure the abundance of CD81 protein specifically associated with exosomes in biological fluids including urine, saliva, cell culture medium, cerebral spinal fluid and blood plasma.

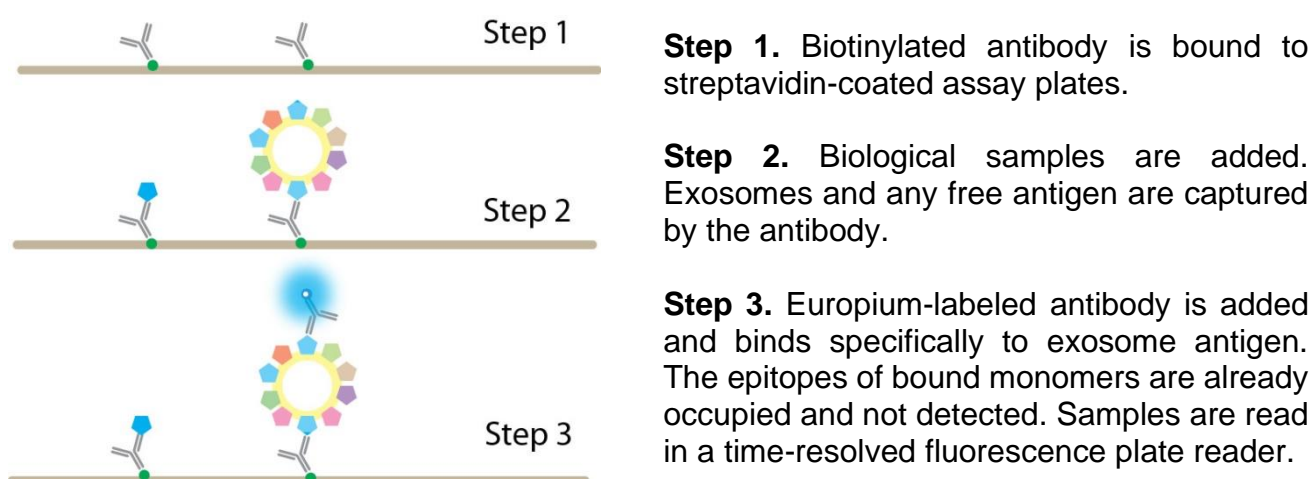


Figure 1. Protocol overview

General note

The technical positive control provided with the kit can be used to verify that adequate technical and handling procedures have been conducted and a signal can be detected at the end of the protocol. The provided technical positive control cannot be used to generate a quantification curve for directly quantifying exosomes. Extracellular vesicles of different origin have different CD81 prevalence on their surface.

If absolute quantification is required, an internal quantification curve can be obtained for each population of exosomes studied by correlating CD81 ExoLISA signal intensity with Nanoparticle Tracking Analysis (NTA) measurements at different concentrations of purified exosomes.

Reagent preparation

- **Preparation of the wash buffer:** Dilute the wash buffer concentrate 25x in Milli-Q water (20 ml concentrate in 500 ml Milli-Q water). The diluted solution may then be stored at room temperature.
- **Preparation of the biotin CD81 in assay buffer:** [Dilutions should be made on the day the assay is performed]. Briefly centrifuge the vial containing the antibody to gather all the liquid at the base of the vial. Prepare a ~2 ng/μl working solution of biotinylated Ab, by diluting 22 μl (22 μg) in 11 ml of assay buffer.
- **Preparation of the europium CD81 in assay buffer:** [Dilutions should be made on the day the assay is performed]. Briefly centrifuge the vial containing the antibody to gather all the liquid at the base of the vial. Prepare a ~0.25 ng/μl working solution of Eu-labeled CD81 mAb, by diluting 40 μl (2.75 μg) in 11 ml of assay buffer.

Technical positive control preparation

1. 25 μg of exosomes have been purified from cell culture and their total protein concentration determined by BCA assay. The 25 μg of exosomes are lyophilized and need to be reconstituted with 250 μl of Milli-Q water to generate a final concentration of **100 μg/ml**.
2. Prepare 7 microcentrifuge tubes, each with 125 μl of PBS.
3. Prepare 8 serial dilutions as shown in the table below. Ensure that samples are properly mixed in each tube.
We recommend beginning with a concentration of 100 μg/ml, as this will ensure that the most diluted sample is within the limit of detection of the assay.

Table 1. Positive control preparation

Technical positive control #	Exosomes	Diluted in PBS (μl)	Exosome concentration (125 μl solution)
1	25 μg (resuspended in 250 μL MQ water)	0	100 $\mu\text{g}/\text{ml}$
2	125 μl of Tube #1	125	50 $\mu\text{g}/\text{ml}$
3	125 μl of Tube #2	125	25 $\mu\text{g}/\text{ml}$
4	125 μl of Tube #3	125	12.5 $\mu\text{g}/\text{ml}$
5	125 μl of Tube #4	125	6.3 $\mu\text{g}/\text{ml}$
6	125 μl of Tube #5	125	3.1 $\mu\text{g}/\text{ml}$
7	125 μl of Tube #6	125	1.6 $\mu\text{g}/\text{ml}$
8	125 μl of Tube #7	125	0.8 $\mu\text{g}/\text{ml}$

See example assay reading on page 9.

Protocol

Coat the wells with CD81 capture antibody

1. Add 100 μl of the freshly prepared dilute solution (2 $\text{ng}/\mu\text{l}$) of biotin-CD81 antibody (prepared as described above) to each well.
2. Incubate the plate for 1 hour at room temperature on a plate shaker at 750 RPM.
3. Wash the plate using an automatic plate washer. Wash each well three times using 250 μl wash buffer for each cycle.
4. Remove the remaining wash buffer.

Add the sample

5. Clear cells and cellular debris from test samples by centrifuging at 3,000 x g for 20 minutes.
6. Transfer 100 µl of the test sample supernatant to each well. Use 100 µl of PBS instead of sample in order to generate a blank reading.
7. Incubate the plate for 1 hour at room temperature on the plate shaker at 750 RPM.
8. Wash the plate using an automatic plate washer. Wash each well three times using 250 µl wash buffer each time.
9. Carefully aspirate the remaining wash buffer.

Add the europium-labeled CD81 detection antibody

10. Add 100 µl per well of the freshly prepared Eu-labeled CD81 antibody dilution (prepared as described above).
11. Incubate the plate for 1 hour at room temperature on the plate shaker at 750 RPM.
12. Wash the plate using an automatic plate washer. Wash each well three times using 250 µl wash buffer each time.
13. Carefully aspirate the remaining wash buffer.

Signal enhancement and reading

14. Add 100 µl of EFI solution to each well.
15. Incubate the plate for 15 minutes at room temperature on the plate shaker at 750 RPM.
16. Measure fluorescence on a time-resolved fluorescence microplate reader as shown on the table below. **Measurements should be performed in triplicate. Before taking the readings, make sure that the plate reader is set to read the fluorescence at the bottom of the plate. *Please check plate reader manual for optimum value of this variable for reading Europium fluorescence.**

Table 2. Instrument settings

Optic settings	Time resolved fluorescence (TRF)	
	Filters	Excitation: 340 nm Emission: 615 nm
General settings	Number of flashes	200 with flash lamp
	Settling time	0.1 - 0.3 s*
	Integration start/lag time	60 - 200 µs*
	Integration time	400 µs

Technical positive control example

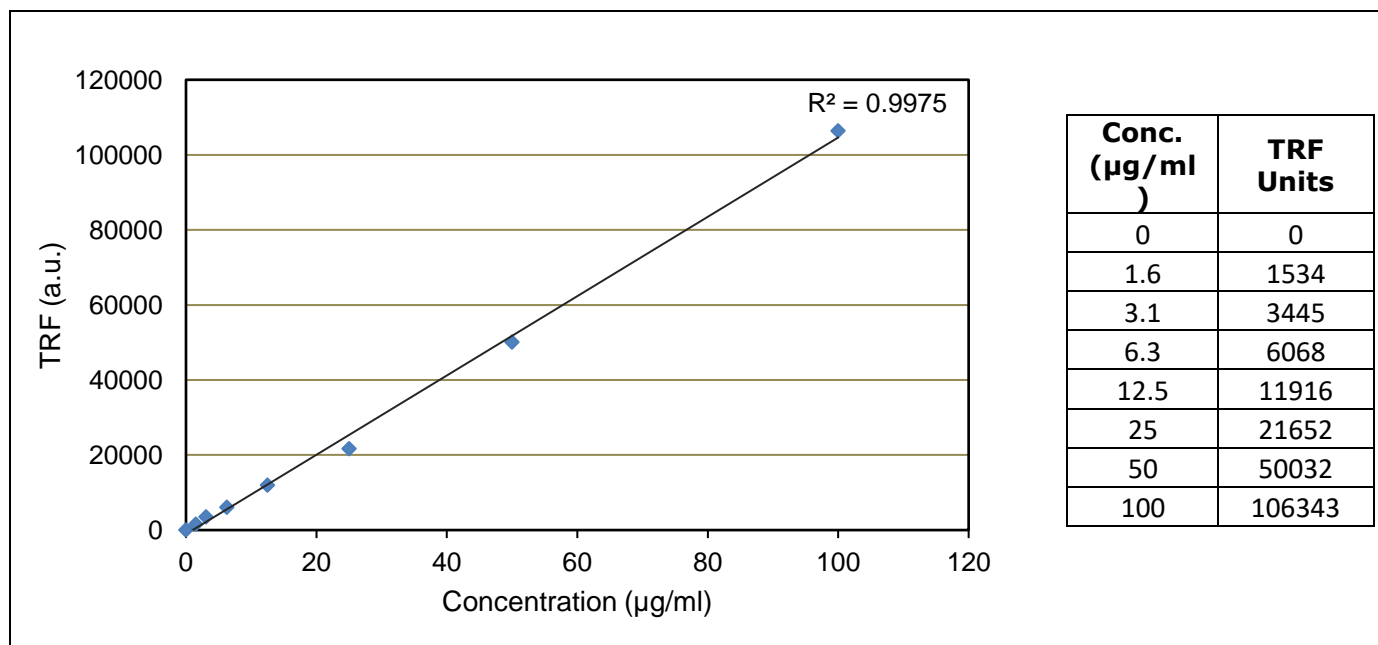


Figure 2. Example readings obtained using the ExoLISA exosome assay with the provided technical control sample (exosomes purified from cell culture). Note the linearity of response over a wide range of concentrations.

Related products

Related products	Product description	Product code
Exosome purification	Exo-spin™ purification kit	EX01, EX02, EX03, EX04, EX05, EX07
Nanoparticle Tracking Analysis (NTA) size profiling service	ZetaView NTA Particle Analysis Service	ZV-1 and ZV-12

Exo-spin™ purification kit

The Exo-spin™ technology combines precipitation and size exclusion chromatography techniques, making it a superior method for exosome separation and concentration, allowing for high specificity and high recovery of exosomes. Exo-spin™ is available in 5 different configurations represented with catalogue codes EX01, EX02, EX03, EX04, EX05, and EX07; specifically designed and optimized for different sample types and downstream applications.

NTA size profiling service

Exosome characterization service for analysis of particle size and particle concentration using the ZetaView instrument from Particle Metrix.

Purchaser Notification

Limited warranty Cell Guidance Systems and/or its affiliate(s) warrant their products as set forth in the Terms of Sale found on the Cell Guidance Systems web site at www.cellgs.com/Pages/Terms_and_Conditions.html

If you have any questions, please contact Cell Guidance Systems.

This product incorporates licensed technologies. The purchase of this product conveys to the purchaser the limited, non-transferable right to use the purchased amount of the product only to perform internal research for the sole benefit of the purchaser. No right to resell this product or any of its components is conveyed. This product is for internal research purposes only and is not for use in commercial services of any kind, including, without limitation, reporting the results of purchaser's activities for a fee or other form of consideration. For information on obtaining additional rights, please contact info@cellgs.com.

CELL GUIDANCE SYSTEMS AND/OR ITS AFFILIATE(S) DISCLAIM ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, OR NON-INFRINGEMENT. TO THE EXTENT ALLOWED BY LAW, IN NO EVENT SHALL CELL GUIDANCE SYSTEMS AND/OR ITS AFFILIATE(S) BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

Cell Guidance Systems' reagents and services enable control, manipulation and monitoring of the cell, both *in vitro* and *in vivo*.

Growth Factors

- Recombinant
- Sustained Release

Exosomes

- Purification
- Detection
- NTA Service

Matrix Proteins

Small Molecules

Cell Counting Reagent

Cytogenetics Analysis



General info@cellgs.com
Technical Enquiries tech@cellgs.com
Quotes quotes@cellgs.com
Orders order@cellgs.com

www.cellgs.com

EUROPE
Cell Guidance Systems Ltd
Maia Building
Babraham Bioscience Campus
Cambridge
CB22 3AT
United Kingdom
T +44 (0) 1223 967316
F +44 (0) 1223 750186

USA
Cell Guidance Systems LLC
Helix Center
1100 Corporate Square Drive
St. Louis
MO 63132
USA
T 760 450 4304
F 314 485 5424