

GFM40 Recombinant Mouse MCP-3 / CCL7

Description

Monocyte Chemotactic Protein 3 (MCP-3), also known as CCL7, is produced by macrophages and tumor cell lines. MCP-3 signals through the G protein-coupled receptors CCR1, CCR2, and CCR3. MCP-3 chemoattracts monocytes and regulates macrophage function during inflammation and metastasis.

Length	74 aa
Molecular Weight	8.5 kDa
Source	E. coli
Accession Number	Q03366
Purity	≥95% determined by reducing and non-reducing SDS-PAGE


Specifications

Alternative Names	Monocyte Chemotactic Protein 1, CCL2, JE, MCAF
Biological Activity	Mouse MCP-3 is fully biologically active when compared to standard. The activity is determined by the ability to induce chemotaxis of THP-1 cells.
Endotoxin Level	≤1.00 EU/μg as measured by kinetic LAL
Formulation	Lyophilized from a sterile (0.2 micron) filtered aqueous solution containing 0.1% Trifluoroacetic Acid (TFA)
AA Sequence	QPDGPNASTC CYVKKQKIPK RNLKSYRRIT SSRCPWEAVI FKTKKGMEVC AEAHQKWVEE AIAYLDMKTP TPKP

Preparation and Storage

Reconstitution	Centrifuge vial before opening. When reconstituting the product, gently pipet and wash down the sides of the vial to ensure full recovery of the protein into solution. It is recommended to reconstitute the lyophilized product with sterile water at 0.1 mg/ml, which can be further diluted into other aqueous solutions.
Stability and Storage	12 months from date of receipt when stored at -20°C to -80°C as supplied. 1 month when stored at 4°C after reconstituting as directed. 3 months when stored at -20°C to -80°C after reconstituting as directed.

Data

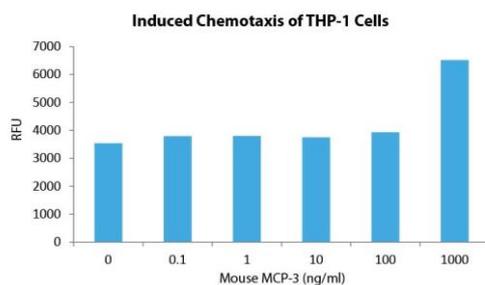

Induced chemotaxis of THP-1 cells assay for Mouse MCP-3. Cells that migrated were counted using a luminescent substrate. Migration over basal levels was reported in response to Mouse MCP-3 starting at 1000 ng/ml.